

aneavTM

Ceramic Bracket System



Trusted aesthetics.


HIGHLAND
metalsTM



Aneav™, the latest addition to Highland Metals™, is a clear bracket designed to deliver durability, comfort and clinical efficiency.

Manufactured using advanced Ceramic Injection Molding (CIM) technology and crafted from high-quality polycrystalline aluminum oxide, Aneav™ provides excellent transparency and performance throughout treatment.


Comfort and Aesthetics

- Low profile design featuring **smooth surfaces** and **rounded edges** for patient comfort.

- Standard ID dot for **bracket identification**.

- **Stain-resistant material** for a pristine, long-lasting look.


- **Placement Guides** assist in alignment with the tooth's long axis.

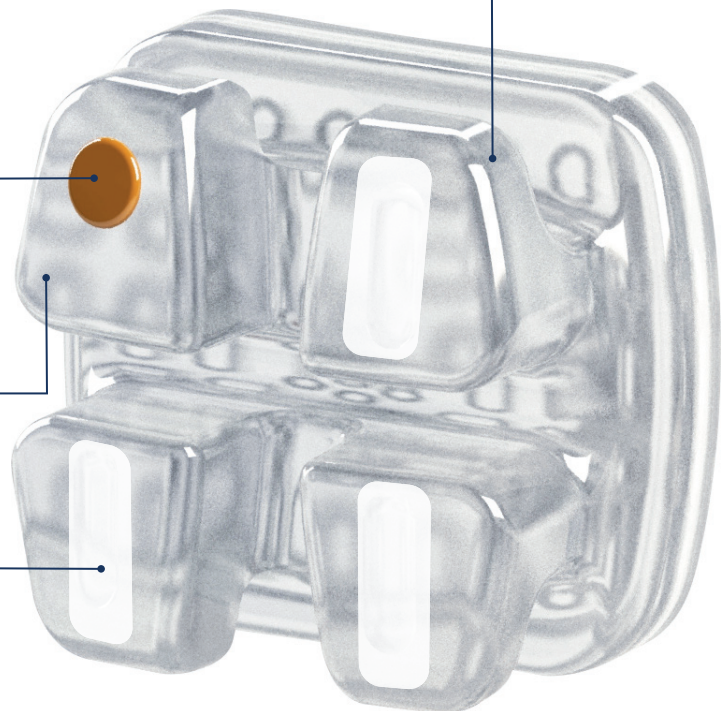

Expansive under tie wing area for **rapid and secure ligation**.

Bond Strength and Fit

- Exclusive lasered texture for **reliable bond strength**.
- **Anatomical mesiodistal and occlusogingival base** curvatures tailored for a precise fit.
- Retention corner in the base with integrated retention groove facilitates **removal of excess adhesive**.

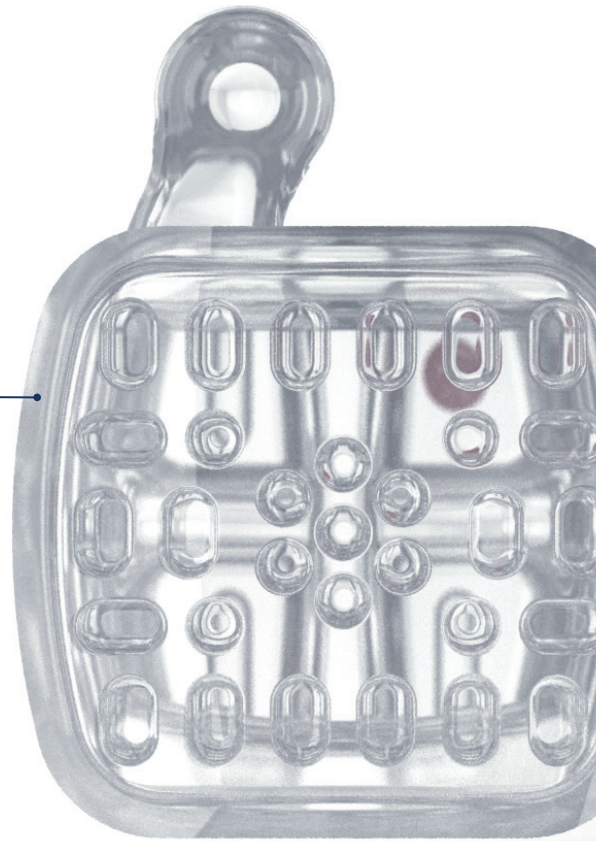

Precision and Performance

- Cuspid and premolar hooks, engineered for durability, **facilitate elastic attachment**.
- Slot with **minimal contact** surfaces **reduces friction**.
- **Built-in torque in base** design simplifies wire mechanics to support a true, straight wire system.





Available in .018 & .022 Roth or MBT¹ prescriptions*

Aneav™ Ceramic Bracket Kits | 20 brackets per kit – upper & lower 5-5

Patient Kits - MBT ¹ Rx	.018	.022
Upper/Lower 5x5 Hook on 3, 4, & 5	135K3451801	135K3452201
Upper/Lower 5x5, 0° Torque Canines (3s), Hook on 3, 4, & 5	135K03451801	135K03452201

Patient Kits - Roth Rx	.018	.022
Upper/Lower 5x5 Hook on 3, 4, & 5	134K3451801	134K3452201

Individual brackets – sold by tooth, in a pack of 5

MBT¹ Rx

	TOOTH	HOOK	TORQUE	ANGLE	ROTATION	ID	RIGHT/LEFT	.018	.022
MAXILLARY	Central		17°	4°		●	R L	135U18RV 135U18LV	135U12RV 135U12LV
	Lateral		10°	8°		●	R L	135U28RV 135U28LV	135U22RV 135U22LV
	Cuspid	Hk	-7°	8°		●	R L	135U38RKV 135U38LKV	135U32RKV 135U32LKV
	Cuspid	Hk	0°	8°		●	R L	135U38RKOV 135U38LKOV	135U32RKOV 135U32LKOV
	Bicuspid	Hk	-7°	0°		●	R L	135U48RKV 135U48LKV	135U42RKV 135U42LKV

	TOOTH	HOOK	TORQUE	ANGLE	ROTATION	ID	RIGHT/LEFT	.018	.022
MANDIBULAR	Anterior		-6°	0°		●	U	135L18UV	135L12UV
	Cuspid	Hk	-6°	3°		●	R L	135L38RKV 135L38LKV	135L32RKV 135L32LKV
	Cuspid	Hk	0°	3°		●	R L	135L38RKOV 135L38LKOV	135L32RKOV 135L32LKOV
	1st Bicuspid	Hk	-12°	2°		●	R L	135L48RKV 135L48LKV	135L42RKV 135L42LKV
	2nd Bicuspid	Hk	-17°	2°		●	R L	135L58RKV 135L58LKV	135L52RKV 135L52LKV

Roth Rx

	TOOTH	HOOK	TORQUE	ANGLE	ROTATION	ID	RIGHT/LEFT	.018	.022
MAXILLARY	Central		12°	5°		●	R L	134U18RV 134U18LV	134U12RV 134U12LV
	Lateral		8°	9°		●	R L	134U28RV 134U28LV	134U22RV 134U22LV
	Cuspid	Hk	-2°	11°	4°	●	R L	134U38RKV 134U38LKV	134U32RKV 134U32LKV
	Bicuspid	Hk	-7°	0°		●	R L	134U48RKV 134U48LKV	134U42RKV 134U42LKV

	TOOTH	HOOK	TORQUE	ANGLE	ROTATION	ID	RIGHT/LEFT	.018	.022
MANDIBULAR	Anterior		0°	0°		●	U	134L18UV	134L12UV
	Cuspid	Hk	-11°	5°	2°	●	R L	134L38RKV 134L38LKV	134L32RKV 134L32LKV
	1st Bicuspid	Hk	-17°	0°		●	R L	134L48RKV 134L48LKV	134L42RKV 134L42LKV
	2nd Bicuspid	Hk	-22°	0°		●	R L	134L58RKV 134L58LKV	134L52RKV 134L52LKV

* HMI™ prescriptions are not implied to be an exact version of any other system, nor do we claim any endorsement by the doctor.

¹ MBT is a trademark of 3M Unitek Corporation.



highlandmetals.com/aneav

Highland Metals
40 Linville Way - Franklin, IN 46131 USA
+1 408-271-2955 (International)

